



Lightweight!

| Weight Innovation | Dual Safeguard | Artificial Intelligence |

EzRay **Air**

“Weight Innovation!”



It all starts with **Weight Innovation**

A light portable intraoral X-ray creates a stable condition for positioning the device. The ease of device positioning produces clear images without any motion artifact.

Handling a lightweight device is very convenient as it reduces the time taken to process a full mouth series. It also produces faster workflow and less stress for your staff.



Incisor



Molar (65 kVp, 2.5 mA, 0.2 sec)



Double scatter shield in One



Internal Shielding

The internal radiation shielding is perfectly designed to protect the operator from radiation leakage .

External Backscatter Shielding

Scatter radiation results from the beam interacting with the surface of the patient, causing radiation to bounce off as residuals scatter in different directions.

The backscatter shield significantly reduces the operator's radiation dose.

“How Clear Image!”




Save your time with **Smart Positioning**

The smart positioning makes preparation simple.

There are default exposure settings for the angle of each tooth.

Operators only need to set the guideline following an occlusal plane or a frankfurt line as the angle and the tooth are automatically detected. The exposure time is also automatically adjusted with the aid of smart positioning.

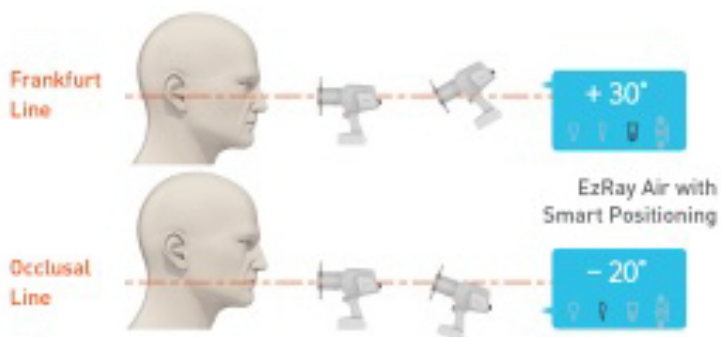

Full Mouth Series within only 2 minutes



The fastest and easiest way to get a full mouth series



One dial for all functions


With EzRay Air, there is no need to remember and press numerous buttons for each setting.

The Smart Dial has all the settings in one dial, reducing the preparation time and creating a simpler and faster work flow.

The Smart Positioning function also provides optional tooth mode settings and scan time settings.



Dimensions [Unit: mm]

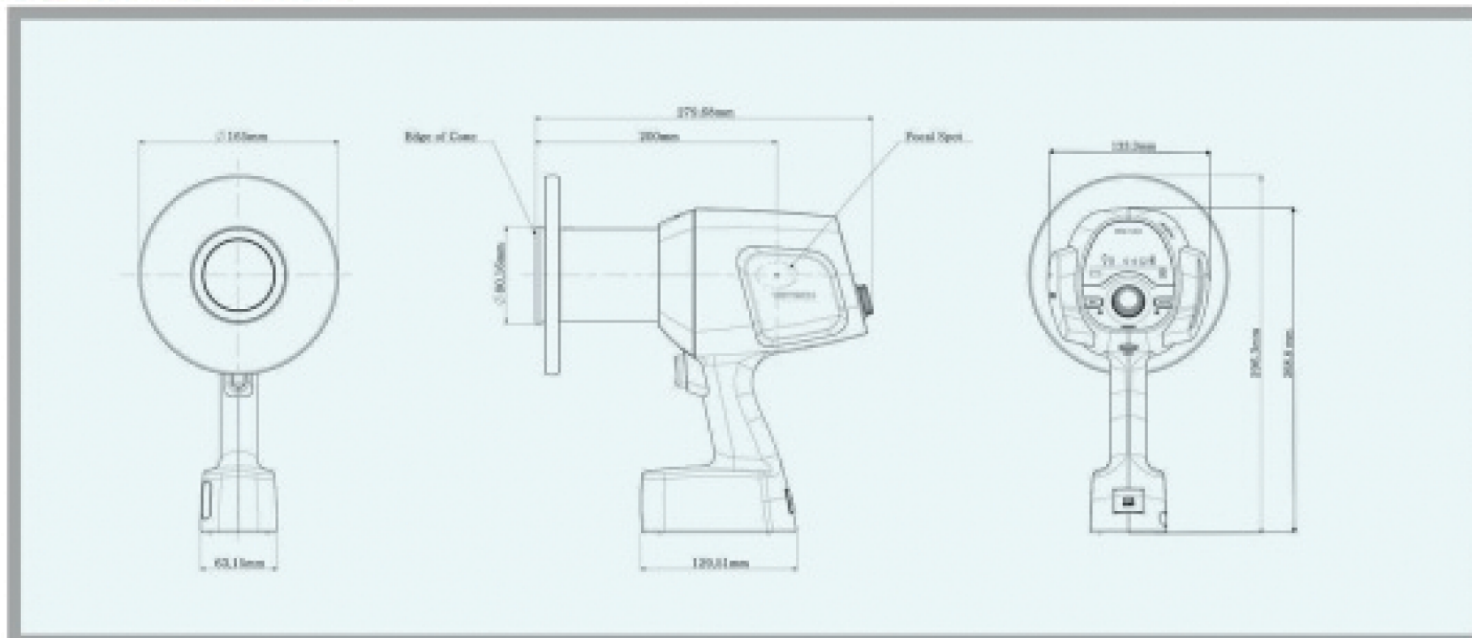

Specifications [EzRay Air : VEX-P300]

Focal Spot	0.4 mm [IEC 60336]
Tube Voltage [kV]	65 kV
Tube Current [mA]	2.5 mA
Exposure Time Range	0.05 – 0.5 sec
Total Filtration	Min. 1.5 mm Al
Source to Skin Distance	200 mm
X-ray Field	Default: 60 mm Round, 30 x 40 mm Rectangular / Optional: 20 x 30 mm Rectangular
Maximum duty cycle	1 : 30
Power Input	22.2 V
Weight	1.4 kg [2.2 lbs.] * without backscatter shield

* The specifications are subject to change without prior notice.